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## User's Manual Fibre-Optic Laryngoscopes

Dear Customer, we thank you for choosing a KaWe product. Our products are known for their high quality and long life. This KaWe product meets EC standards 93/42/EEC (standards for medical products). All of our fibre-optic blades are made of stainless steel DIN 1.4301/1.4305 in accordance with ISO 7376-3 standards.

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### Purpose:

Used for direct inspection of the larynx. These laryngoscopes are mainly used in cases of medical emergency for oral intubation of accident victims or for intubation before inducing anaesthesia during operations. They are used in connection with an endotracheal tube in order to remove acute respiratory disturbances and before operations in hospitals in order to ensure a clear respiratory tract. This instrument consists of 3 principal parts: the handle, the light source and the blade.

### Unsuited Use/Contraindication of the Products:

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### Risks / Warning of possible dangers:

- Danger of injury of the sensible mucosa and the epiglottis during the intubation.
- During intubation, ensure that the lips, the teeth, the tongue, the larynx and the vocal cords are not in danger of being harmed. In most cases, injuries to the lips and the tongue are caused by confinement, whereas injuries to the teeth are mostly caused by use of the laryngoscope.
- Always use the laryngoscope blades must be disinfected and sterilized before further use (autoclave treatment, however, thermolabile media - beware towards shortened service life with bright guide due to thermal and mechanical stress).

After use, please clean the handles with an alkaline disinfectant solution in order to avoid possible infection.

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#### User's Manual Fibre-Optic Laryngoscopes

Dear Customer, we thank you for choosing a KaWe product. Our products are known for their high quality and long life. This KaWe-product meets EC standards 93/42/EEC (standards for medical products). All of our fibre-optic blades are made of stainless steel DIN 1.4301/1.4305 in accordance with ISO 7376-3 standards.

Please read these instructions thoroughly and carefully before attempting to use this product and heed the given instructions.

Utilization: The laryngoscopes are only to be used by authorized personnel.

#### Purpose:

Used for direct inspection of the larynx. These laryngoscopes are mainly used in cases of medical emergency for oral intubation of accident victims or for intubation before inducing anaesthesia during operations. They are used in connection with an endotracheal tube in order to remove acute respiration disturbances and before operations in hospitals in order to ensure a clear respiratory tract. This instrument consists of 3 principal parts: the handle, the light source and the blade.

#### Unsuited Use/Contraindication of the Products:

Any use other than that described here is not in accordance with the intended use of the instrument. The manufacturer is not liable for any resulting damages. The user alone bears the risk.

#### Risks / Warning of possible dangers:

- Danger of injury of the sensitive mucosa and the epiglottis during the intubation.
- During intubation, emphysema of the lungs, the trachea, the larynx and the vocal cords are not in danger of being damaged. In most cases, injuries to the lips and the tongue are caused by confinement whereas injuries to the teeth are mostly caused by use of the laryngoscope.
- After use, the laryngoscope blades must be disinfected and best sterilized before further use (autoclave treatment, however, this must be preceded by a thorough cleaning of the light guide due to thermal and mechanical stress).

- After use, please clean the handles with an alcoholic disinfectant solution in order to avoid possible infection.

#### Maintenance instructions:

##### • Replacing the fibre-optic light guide

Please note that the "MEGALIGHT" blades with integrated fibre-optics cannot be replaced!

##### • Replacing the lamp in the fibre-optic handle

Hint: The satisfactory operation of the laryngoscopes can only be guaranteed when suitable xenon bulbs (REF. 12.75141.003) are used. Remove the head from the handle sleeve. Unscrew the integrated inner enclosure from the inside part of the head! Next, remove the protective cap of the lamp and unscrew and remove the lamp. If necessary, clean the glass bulb of the new lamp with alcohol. The glass bulb must be properly clean and free from fingerprints (free of grease). The new lamp must be screwed completely into the socket.


\* Saw cutters for easy replacement of the light source.

#### Further information, maintenance, storage:

- Familiarize yourself thoroughly with the assembly of the laryngoscopes when they are delivered.
- In the case of frequent intubations, the electric bulb in the fibre-optic handle may become loose and vibrations of the light source may result. Therefore, always ensure that the electric bulb has been inserted correctly!
- After use, disassemble and clean the fibre-optic blades with exchangeable light guides. REF. 03.22013.602 – 03.22013.642 it is recommended that the screw on the light guide be adjusted at regular intervals.

Used and stored properly, this instrument will serve you many years. Pay attention to the following instructions:

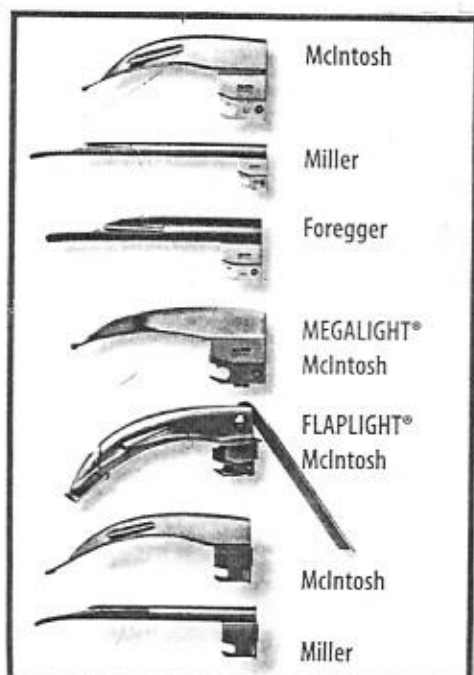medical products made of stainless steel and glass fibres. In addition, the cleaning agents used for economy series blades must be suitable for medical synthetic products. Please ensure that a neutral or alkaline cleaning agent is used. The manufacturer's instructions for the used cleaning agents and of the cleaning machine must be strictly observed.

#### Cold sterilization / Disinfection:

The blades can be disinfected in a 2% glutaraldehyde solution (CIDEK, SPORICIDIN). When cleaning with chemicals, pay careful attention to the instructions (user's manual) given by the individual manufacturers with regard to cleaning time and concentration!

**Attention: Do not clean the blades in chlorine, betadine or peroxide solutions!**

**Gas sterilization:** Chemical / thermal disinfections (ethylene oxide) up to 65°C may be used.

#### Steam sterilization:

Steam sterilization in autoclaves up to 134°C is possible; however, this may shorten the life time of the light guide. We recommend not exposing the light guides to these temperatures as the life time of the fibre-optics will be significantly reduced!

**Ultrasonic cleaning is not allowed!! Flash autoclave treatment and hot air sterilization will damage the blades!**

#### Preparation of the battery-charging handles:

All parts of the handle and of the bulb cap can be wiped off with a cloth moistened with alcohol and with a mild cleaning agent. In addition, the handle sleeve can be treated in the autoclave. Never bring handles with batteries and/or accumulators into contact with water or solutions because these parts have electrical contacts.

#### Attention:

3.5V battery handles that can be recharged in the socket have a separate user's manual (REF 03.21000.812-822).

Should you need further information, please contact your dealer!

Manufacturer: KaWe

#### Warning:

Please pay attention to the manufacturer's instructions! Never bring the laryngoscope into contact with liquid and make sure that no liquid penetrates the handle! Never lay the blades into chlorine, betadine or peroxide solutions. Sterilization in hot air devices is not allowed. In order to prevent crevice corrosion or rust, never bring the laryngoscopes into contact with salt solution. This may result in pitting corrosion.

#### Limitation of reuse preparation:

Steam sterilization up to 134°C is possible, but may have an effect on the life time of the fibre-optics.

#### Instructions:

##### Workstation:

It is recommended that the blades be immediately cleaned with running water or a soap solution after their use.

##### Storage and transport:

In order to avoid dirt accumulation and dust, always store the laryngoscope in its packaging.

##### Cleaning preparations:

Remove the blade from the handle and follow the given cleaning instructions.

##### Manual cleaning:

After removing the batteries/accumulators and the bulb cap, the handles can be wiped off with a fabric moistened with alcohol. Clean the blades with running water or a soap solution.

##### Automatic cleaning:

Automatic cleaning of the blades with a neutral or alkaline cleaning agent. For automatic cleaning, please always use special preparation systems with the intended use of this product.

##### Disinfection:

Disinfect the blades with a 2% glutaraldehyde solution (CIDEK, SPORICIDIN). If a chemical solution is used, the manufacturer's instructions (user's manual) must be followed. After disinfection, clean blades with clear and running water in order to prevent any residues from drying onto the material.

##### Maintenance:

Defective parts must be removed and/or replaced. Defective parts must be exchanged. Check for the tight fit of the bulb. We recommend tightening the light guide screw prior to each use. Defective light guides must be replaced. Prior to each use, ensure that the units are completely serviceable.

##### Check and performance test:

Prior to each use, check the performance of the lamp, the batteries and/or accumulators. Defective units and units that are to be repaired should be stored separately in order to prevent their reuse. In case that a unit is returned to KaWe, proof of a previous disinfection or sterilization must be provided.

**Packaging:** Individual packaging in bags or cartons.

##### Sterilization:

Chemical/thermal gas sterilization of blades up to 65°C is allowed as is steam sterilization in autoclaves up to 134°C (steam quality in accordance with EN 285). The handle sleeves can be sterilized in the autoclave up to 134°C (without batteries and/or accumulators and bulb caps). Please pay attention to the instructions of the individual manufacturers of the sterilizers. Ensure that the filling weight of the autoclave is not exceeded. The blades / handle sleeves must be dry when they are treated in the autoclave, or in the last reuse cycle prior to their treatment in the autoclave, they must be cleaned with distilled water.

**Storage:** Store the dry units in a clean room at ambient temperature.

##### Additional information:

Please pay careful attention to the instructions in the user's manual. In case of doubt, please contact the manufacturer or the dealer.

##### Contacting the manufacturer:

Address or phone number of the dealer or at the number (+49-7141-68188-0)

\*The following endorsed cleaning agents have been used on this product: NEODISHER® MEDICLEAN, Mucapur® CD and Mucapur® XL

##### Guarantee:

When used under normal circumstances and with attention to these instructions (the user's manual) we guarantee this product for two years after the date of purchase (with exception to wear parts such as batteries / accumulators and lamps). Should you need further information or should your instrument require repair, please contact your dealer!

Information about spare and accessory parts can be obtained in our KaWe price list or in our Internet-Shop at [www.kawemed.de](http://www.kawemed.de)

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